



SWISS MADE
SWISS TECHNOLOGY

SERIES 520
ADVANCED/STANDARD MODELS



PRECISA DNA: ZERO COMPROMISE

Laboratory and industrial technicians are aware of the significance of the results they provide. That is why they are unwilling to accept any compromise in their field of work. This is precisely why leading specialists around the world place their trust in Precisa and have done for decades.


PRECISA DNA: SWISS PRECISION

Precisa products are prime examples of Swiss perfection and reliability. They are tried and tested in both demanding laboratory environments and industrial applications, and are subject to rigorous quality controls throughout the entire manufacturing process. Precisa invests continuously in developing new technologies and employs a highly qualified team. The result is a state-of-the-art range of precision instruments for the most exacting demands.

PRECISA DNA: SERIES 520

Swiss idea, Swiss development, Swiss production. With over 50 years knowledge and experience, Precisa Gravimetrics AG has developed one of the most important innovations in weighing technology in recent years.

The Precisa PHASTbloc™ weighing cell allows you to benefit from this advanced technology and provides a technically perfect, easy-to-use, attractive, economical and sustainable balance. Discover the innovations of the 520 series. We look forward to inspiring you - as many users have been already.

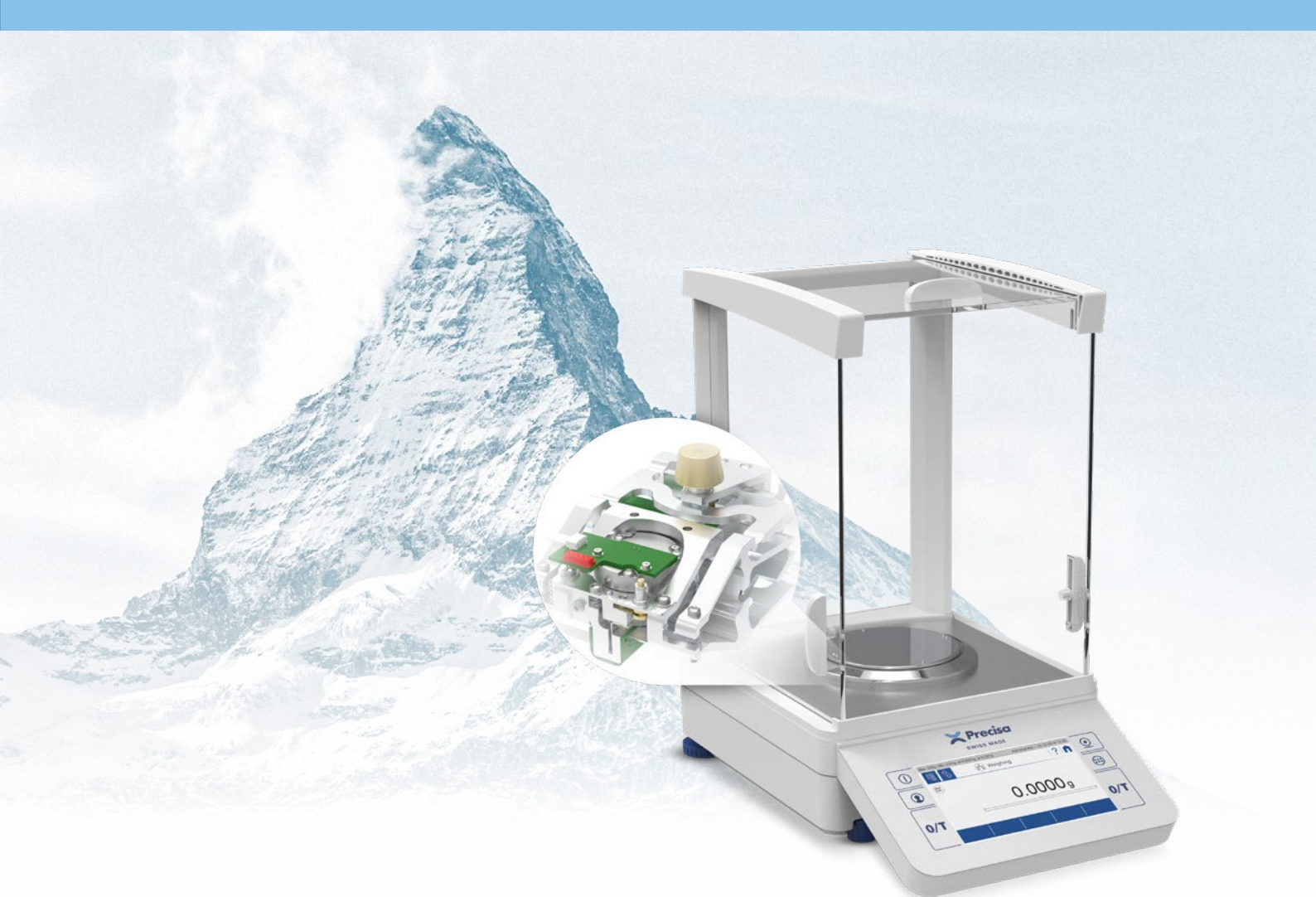


PHASTbloc™


PHASTbloc™ (Precisa Hybrid Advanced Sensor Technology) - Precisa's new and unique weighing technology combines the benefits of current weighing cells with monolithic design. The ultimate in Swiss engineering from Precisa has created a world first weighing technology, representing 50 years of knowledge and expertise in developing weighing technologies and innovations. The sustainable hybrid approach in the PHASTbloc™ unites the benefits of a conventional weighing cell and monolithic design. This enables the efficient maintenance and repair of a weighing cell instead of replacing it - sustainability in practice. Advanced materials, engineering and unique manufacturing methods lift quality and performance to new heights.

Using a range of materials such as standard grade and high tech aerospace aluminium, maximises reparability and performance, offering affordability and an environmentally friendly product. All parts of the PHASTbloc™ are manufactured with the highest Swiss quality and precision. The division of the weighing cell into functional parts not only allows an optimal choice of material and manufacturing process for each section, but also opens up the possibility of replacing parts. This allows the PHASTbloc™ to be maintained cost-effectively and ensures a long life.

All this makes the PHASTbloc™ weighing cell the most ecological and economical of its kind without compromising on performance.

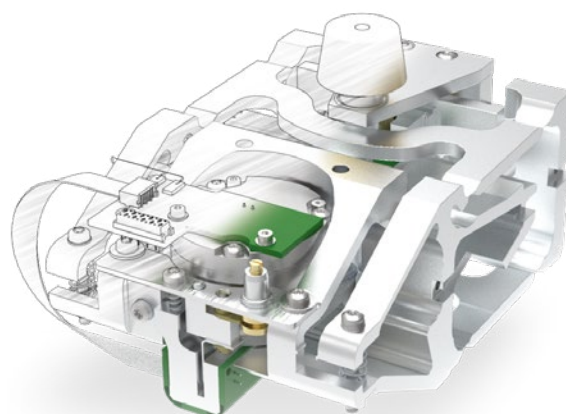

520 PT FEATURES

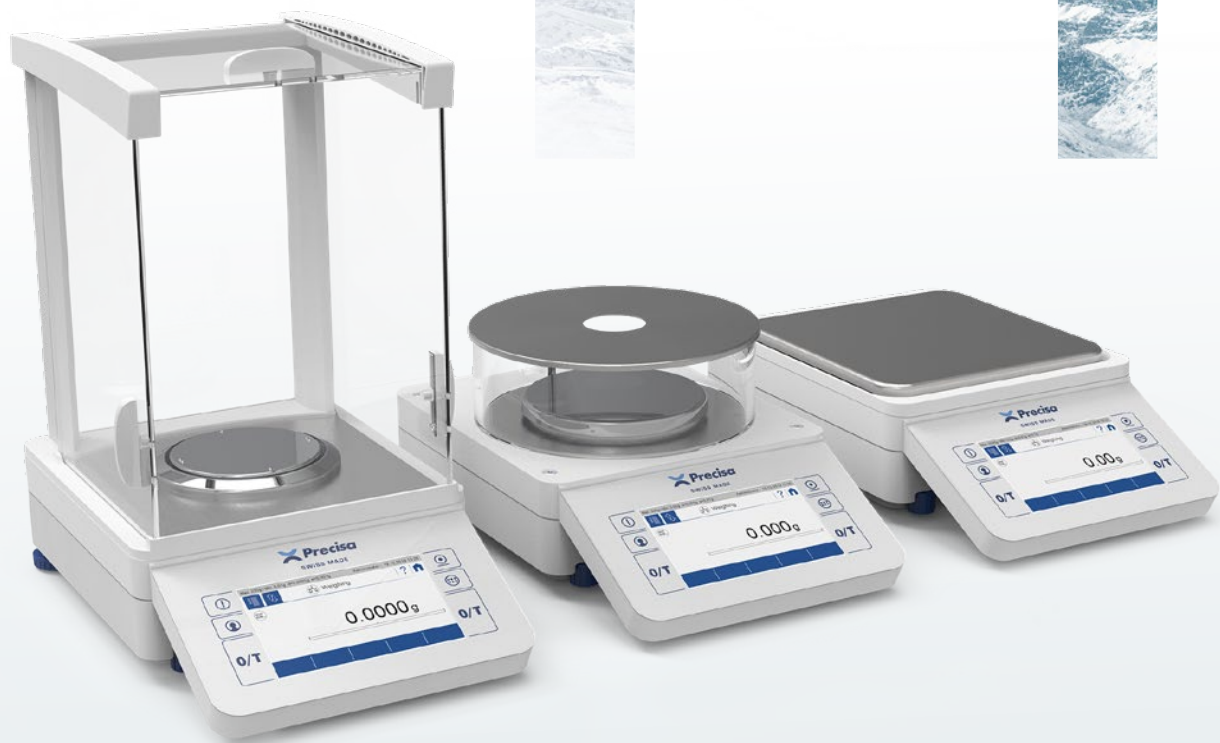

520 PT: THE SMART BALANCE

With a 6.4" touch control screen, this is the largest in its class! In working environments such as laboratories, research and development, where flexibility and individual user solutions are required, the 520 PT reveals its full potential. With the advantages of the touch screen in terms of visualization, intuitive operation and customization options, the 520 PT offers the highest degree of flexibility, safety and user-friendliness. Six large fixed keys for the most common operating functions - on/off switch, tare, printing, user and application selection - optimize operation without reducing the display area.

With a complete range of applications such as counting, density, animal weighing, percentage, check, ... and the RS232, USB, plus optional Ethernet, WiFi and Bluetooth interfaces, there are no limits to the 520 PT. The 520 PT with its touch screen is the right choice for maximum user-friendliness and flexibility.

Features:



Optional / On Demand:

Applications:



Optional / On Demand:

See Features/Applications for an explanation of terms.

520 PT FEATURES



Touchscreen



Applications



Integrated Help Functionality

PRECISA SERIES 520 PT

Touchscreen

Large 6.4" touch-sensitive user interface with an elegant and bright 5" colour display.

As with mobile devices, the balance's user interface is intuitive and easy to use. Precisa's user-friendly operating concept makes navigation and configuration much easier thanks to the colour scheme.

The desired application is just a fingertip away. The colour logic displays all touchable surfaces in blue. The "Home" button, which is always available, guides you back to the active application.

- 6.4" touch-sensitive user interface
- 5" anti-reflective display
- Shock and scratch resistant
- Suitable to operate with laboratory gloves
- UV resistant
- Frameless for easy and thorough cleaning.

6 Fixed Keys

Six quick-select keys with the most important functions.

To the right and left of the display there are three fixed position touch-sensitive keys. This gives you access at all times to the most important functions: switch on/off, print, application and user screen as well as taring or zeroing.

Context Help

We simplify your life!

Simply touch the help icon and the context help appears on screen. You will find all relevant information, even if you are in an application, the context menu, the user screen or anywhere else. Help is available at the touch of a button! No more tedious and time-consuming browsing through user manuals.

520 PT FEATURES



Customizable User Interface



Manage Users and Rights-Access



Electronic Levelling Indicator

PRECISA SERIES 520 PT

Customize Your Screen

Individual and efficient.

The display of the 520 PT series can be fully adapted to suit your requirements. You decide what information is displayed in the notification area.

Six soft buttons can be freely configured in the action bar. Assign the most important functions to them or create a folder with up to four commands. Complex functions and processes are simplified in this way so that you can concentrate on your main tasks.

Nothing stands in the way of simple, efficient and safe work.

User Management

Individual and secure.

With the built-in, extensive user management, your balance can be optimally adjusted for up to 25 users.

Simply select your profile at the start of work and the balance loads your personal settings.

The balance's individual configuration is saved in your user profile: whether with full access rights, restricted to a single application or simply with the basic operating settings for right-handed or left-handed users. Each profile can be protected with a unique password.

Electronic Level Indicator

Improved precision.

With the built-in electronic levelling, you can be assured that your balance is level at all times. If the balance is no longer level, a message appears on the display. The balance can also be configured so that it can only be operated when level.

To reset the position of the balance correctly, a levelling guide is displayed on the screen. It shows which direction the feet need to be moved.

520 PT FEATURES



Communication and Interfaces

Unlimited connectivity.

The 520 PT Series meets your communication needs. Whether you already have a LIMS in a paperless laboratory or not, we offer the balance with a range of standard interfaces:

- USB type A (host)
- USB type B (device)
- RS232
- 10/100M Ethernet (optional)
- Wi-Fi (optional)

The integration of your balance into the network, and the connection to a PC or printer is very easy.



520 PB FEATURES




520 PB: TIMELESS AND PROVEN IN EVERYDAY PRACTICE

Five high-quality mechanical membrane keys guarantee longevity and a tactile feel with clearly noticeable feedback to the user. The 520 PB is also suitable for laboratories, the pharmaceutical, food, light and consumer goods industries as well as production or for pharmacies and jewellery. Wherever precise and demanding routine work is carried out on a regular basis, reliability and user-friendliness are key components.

State-of-the-art with an extensive range of applications such as counting, density, animal weighing etc. Including statistics, a wide range of languages, RS232 and USB interfaces, all current and future requirements are met. It has never been so easy to choose the right balance.


Features:

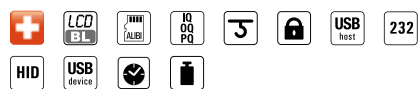

Optional /
On Demand:



Applications:

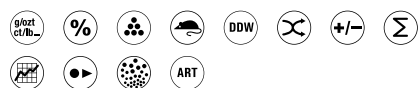

See Features/Applications for an explanation of terms.

520 PB FEATURES

PRECISA SERIES 520 PB

Robust & Durable Design

Designed for almost any environment.

When designing the 520 PB series, attention was paid to robustness and durability. The housing materials, keyboard and display were chosen to withstand rough environments and meet heavy demands.

Large Alpha Numeric LCD

Large and easy to read.

The bright, high-contrast LC display with backlight and large digits (> 15mm) offers excellent readability.

5 Key Membrane Keypad

Clear and simple.

The five softkeys structure the various application possibilities in a user-friendly manner and considerably improve the operation of the device.

Navigation has been enhanced and simplified for the 520 PB series.

Extensive Range of Applications

Advanced functionality.


The desired application can be easily selected from a large number of specific applications and guarantees maximum functionality. The 520 PB series offers an extremely wide range of applications for almost every possible application.



Robust and Durable Design



5 Key Membrane Keypad

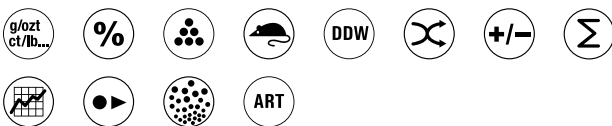


Extensive Range of Applications


520 COMMON FEATURES



Sturdy Housing

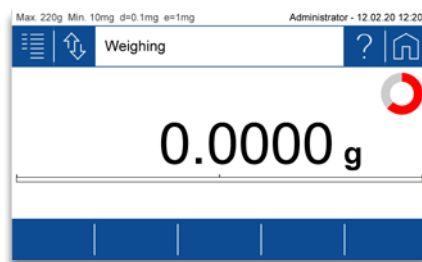

Real-time Stability Indicator

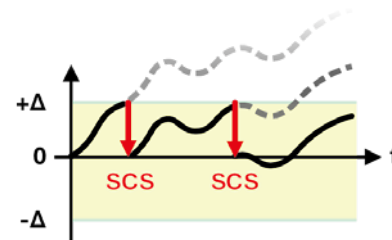

Self-calibration System

Sturdy Housing

Robust and durable.

Precisa's instrument housings have always been characterized by their robustness and longevity. This has not changed with the current Series 520.

Both the new PHASTbloc™ weighing cell and the entire balance have been designed for optimal material choice. For this reason, the base is a robust die-cast aluminium housing and the upper part is made of fibreglass-reinforced and durable thermoplastic compound.

Stability Indicator

Real-time stability indicator with stability prediction.


The optical stability indicator indicates measurement faults and provides real-time information on the influences of the environment which may affect the balance.

The more stable the environment, the faster the stability indicator disappears and a reliable value is displayed. The stability indicator allows estimating how long it takes for the balance to display a reliable value.


Self Calibration System (SCS)

Always accurate.


The Self Calibration System (SCS) developed by Precisa enables automatic adjustment / calibration of the balance. The function can be set to be temperature and/or time-controlled, so that the user can always expect maximum measuring accuracy. This option is often used in regulated environments (GLP).

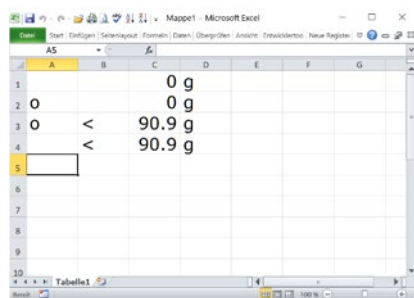

Human Input Device (HID)

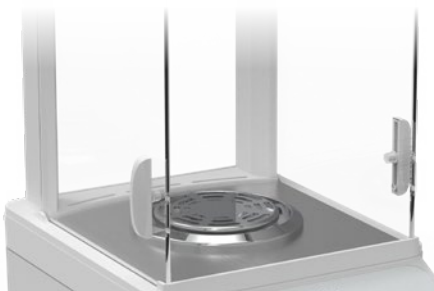

Enhanced Draft Shield

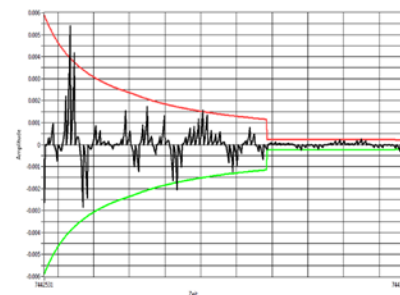

ENC™ Environmental Noise Cancelling System

Human Input Device (HID)

Direct balances – PC communication.

The balance can be connected to a computer via a USB cable and is recognized by the computer as a human input device (HID). This turns your balance into a keyboard and weight values can be sent directly from the balance to the PC without installing any software. No additional authorizations, apps or tools are required. In this way, the data can be written directly into the open program, for example Excel, Word, Notepad, etc.

Enhanced Draft Shield Technology

For the use in safety cabinets.

Employee safety is our highest priority. When working with hazardous samples, weighing must be carried out in appropriate safety enclosures, such as fume hoods or safety cabinets.

The innovative mechanical design of the new draft shield enables precise and reliable weighing, even under laminar airflow conditions.

Other Features

Suitable for a wide range of uses.

- Multilingual
- Improved digital filter to eliminate external effects (ENC™)
- Space saving
- Maximum stability, thanks to four adjustable feet
- Weighing below option
- Optionally IP65 protected

MODELS

PT MODEL

Model	Capacity	Readability	Repeatability (Load >5%)	Min. Weight (1%, k=2, Load >5%)	Linearity	Pan Size (mm)	Class / „e“	Response Time	Item No.
Analytical Balances									Net weight 6.6 kg
PT 120A	120 g	0.1 mg	0.1 mg	20.0 mg*	0.2 mg	ø 90	I / 1 mg**	2 s	520-9221-002***
PT 220A	220 g	0.1 mg	0.1 mg	20.0 mg*	0.2 mg	ø 90	I / 1 mg**	2 s	520-9222-002***
Precision Balances									Net weight 5.6 kg
PT 320M	320 g	1 mg	1 mg	200 mg*	1.5 mg	ø 120	II / 10 mg**	1.5 s	520-9321-002***
PT 620M	620 g	1 mg	1 mg	200 mg*	1.5 mg	ø 120	II / 10 mg**	1.5 s	520-9322-002***
PT 320C	320 g	0.01 g	0.005 g	1 g*	10 mg	186 x 186	-	1.5 s	520-9421-002***
PT 620C	620 g	0.01 g	0.005 g	1 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9422-002***
PT 1200C	1200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9423-002***
PT 2200C	2200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9424-002***
PT 3200C	3200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9425-002***
PT 4200C	4200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9426-002***
PT 6200C	6200g	0.01 g	0.01 g	2 g*	15 mg	186 x 186	-	1.5 s	520-9427-002***
PT 3200D	3200 g	0.1 g	0.05 g	10 g*	0.1 g	186 x 186	II / 0.1 g**	1 s	520-9521-002***
PT 6200D	6200 g	0.1 g	0.05 g	10 g*	0.1 g	186 x 186	II / 1.0 g**	1 s	520-9522-002***
PT 10200D	10200 g	0.1 g	0.1 g	20 g*	0.15 g	186 x 186	II / 1.0 g**	2 s	520-9523-002***
PT 10200G	10200 g	1 g	0.5 g	100 g*	0.5 g	186 x 186	II / 1.0 g**	1.5 s	520-9621-002***

PB MODEL

Model	Capacity	Readability	Repeatability (Load >5%)	Min. Weight (1%, k=2, Load >5%)	Linearity	Pan Size (mm)	Class / „e“	Response Time	Item No.
Semi-Micro Balances									Net weight 6.6 kg
PB 125SM-DR	52 g / 125 g	0.01 mg / 0.1 mg	0.015 mg / 0.05 mg	3.0 mg*	0.06 mg / 0.1 mg	ø 80	–	8 s	520-9112-001
Analytical Balances									Net weight 6.6 kg
PB 120A	120 g	0.1 mg	0.1 mg	20.0 mg*	0.2 mg	ø 90	I / 1 mg**	2 s	520-9211-002***
PB 220A	220 g	0.1 mg	0.1 mg	20.0 mg*	0.2 mg	ø 90	I / 1 mg**	2 s	520-9212-002***
Precision Balances									Net weight 5.6 kg
PB 320M	320 g	1 mg	1 mg	200 mg*	1.5 mg	ø 120	II / 10 mg**	1.5 s	520-9311-002***
PB 620M	620 g	1 mg	1 mg	200 mg*	1.5 mg	ø 120	II / 10 mg**	1.5 s	520-9312-002***
PB 320C	320 g	0.01 g	0.005 g	1 g*	10 mg	186 x 186	–	1.5 s	520-9411-002***
PB 620C	620 g	0.01 g	0.005 g	1 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9412-002***
PB 1200C	1200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9413-002***
PB 2200C	2200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9414-002***
PB 3200C	3200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9415-002***
PB 4200C	4200 g	0.01 g	0.01 g	2 g*	10 mg	186 x 186	II / 0.1 g**	1.5 s	520-9416-002***
PB 6200C	6200 g	0.01 g	0.01 g	2 g*	15 mg	186 x 186	–	1.5 s	520-9417-002***
PB 3200D	3200 g	0.1 g	0.05 g	10 g*	0.1 g	186 x 186	II / 0.1 g**	1 s	520-9511-002***
PB 6200D	6200 g	0.1 g	0.05 g	10 g*	0.1 g	186 x 186	II / 1.0 g**	1 s	520-9512-002***
PB 10200D	10200 g	0.1 g	0.1 g	20 g*	0.15 g	186 x 186	II / 1.0 g**	2 s	520-9513-002***
PB 10200G	10200 g	1 g	0.5 g	100 g*	0.5 g	186 x 186	II / 1.0 g**	1.5 s	520-9611-002***

* USP Min. weight (0.1%, k=2, Load >5%) = Min. weight x 10

** EU/OIML/UKCA type approval on demand, Initially verified with SCS (-001) only

*** -001 for balance with SCS / -002 for balance without SCS

OPTIONS AND ACCESSORIES

Options	PT	PB
EU-Verification according to 2014/31/EU by Precisa (initial verification) incl. alibi memory, available only for balances with EU/OIML type approval and SCS	350-8938	350-8938
520 PB Protection IP65, protection against ingress of dust and low pressure water jets, fitted at factory, not all types available	–	350-8754
520 PT Protection IP50, protection against ingress of dust, fitted at factory, not all types available	350-8748	–
SCS Self calibration System with clock, automatic on temperature change and real-time programmable	520-xxxx-001	520-xxxx-001
Balance with AeroPan, with SCS	520-xxxx-005	520-xxxx-005
Balance with AeroPan, without SCS	520-xxxx-006	520-xxxx-006

Applications	PT	PB
BEST Buoyancy Error Suppression Technology	350-8739-001	–

Accessories	PT	PB
Draft shield PB	350-8735	350-8735
Draft shield PB + Weighing pan M* Includes all necessary parts: 350-8735, Weighing pan M*	350-8735-001	350-8735-001
Draft shield PB short	350-8742	350-8742
Draft shield PB short + Weighing pan M* Includes all necessary parts: 350-8742, Weighing pan M*	350-8742-001	350-8742-001
AeroPan 390 Ø 90 mm for 0.1 mg balances	350-8729	350-8729
AeroPan 520 186 x 186 for C, D, G balances (for retrofitting 1200 g+ capacities, a re-linearization is recommended)	350-8758	350-8758
Remote display free-standing, incl. 1.5 m cable	-	350-8401-001
Remote display with wall bracket, incl. 1.5 m cable	-	350-8401-002
Hook for weighing below the balance	350-8744	350-8744
Animal weighing bowl complete for 0.01 g, 0.1 g and 1 g balances	350-8743	350-8743
Density kit for 0.1mg balances Container size Ø 75 mm, h = 100mm	350-8719	350-8719
Density kit solids only (w/o Glass body and Hook) for 0.1 mg balances Container size Ø 75 mm, h = 100 mm	350-8720	350-8720
Glass body 10 ccm for density determination of liquids	350-7054	350-7054
Downholder for samples with density < 1 g/cm ³	350-7194	350-7194
Single foot switch USB for one function: tare, print, ...	350-8747-011	350-8747-011
Double foot switch USB for two functions: tare and print, tare and ...	350-8747-012	350-8747-012
Protective cover 520 for the whole balance (0.1 mg, 1 mg), set of 5 pieces	350-8745-001	350-8745-001
Protective cover 520 for the whole balance (0.01 g, 0.1 g, 1 g), set of 5 pieces	350-8745-002	350-8745-002
Protective cover 520 for the whole balance (0.01 mg, 0.1 mg), set of 5 pieces	350-8745-003	350-8745-003
Protective cover for terminal, set of 5 pieces	350-8746	350-8746
Carrying case 520	350-8746	350-8746
LAN Ethernet connection 520 PT	350-8939-000	-
WiFi connection 520 PT	350-8941	-
Interface box, LAN ethernet connection	-	350-8940-001
Interface box, WiFi connection	-	350-8940-002

Operating Temperature	Line Voltage, Tolerance	Power Consumption	Frequency
5 .. 40 °C (41°F .. 104°F)	115V / 230V, +15% / -20%	18 VA	50 - 60 Hz

SPECIAL PRODUCTS / ACCESSORIES





Common Accessories	520 PT/PB
Printer LP4024, 230V / 115V	350-8391
Bluetooth connection for Printer LP4024	350-8391-001
Paper roll for Printer LP4024, set of 4 pieces	350-8392
Adhesive paper roll for Printer LP4024, set of 5 pieces	350-8396
Ribbon for Printer LP4024	PN 3953-013
Ribbon for Printer CBM910	350-8367
Paper roll for Printer CBM910	350-8366
Barcodescanner, USB Typ-A	350-8500-001
Data cable DB9m-RJ45, 1.5m, for balance used as reference balance	350-8661
Data cable DB9m-DB9f, 1.5m, (PC/Printer LP 4024, 1:1)	350-8672
Data cable DB9m-DB25m, 1.5m, (Printer CBM910)	350-8673
USB cable A-B, 1.5m, (PC)	PN 3950-125
Interface box, Interface and power splitter	350-8940-000
HID cable, RS232, DB9m-USB-A, 1.5m, (PC), (x: 0=e, 1=d, 2=f)	use PN 3950-125
External USV accu power pack for continuous operation	350-8662
Power splitter cable	350-8399
Anti-theft cable	350-8555
Carat bowl (for diamond weighing)	350-8322

PrecisaBUS Accessories (available on demand only)	520 PT/PB
Signal Lamp, indicator with 3 lamps, green / red / yellow	350-8734
Horn for signal lamp	350-8533
Input / Output module with 6 TTL inputs / 8 relay outputs	on demand

FEATURES & APPLICATIONS

Features Series 520		PT	PB	Applications Series 520		PT	PB
	Swiss Made	✓	✓	 Units g, mg, ozt, lb, tael, etc.		✓	✓
	Color Touch Screen	✓		 Percent weighing		✓	✓
	Backlit LCD		✓	 Piece counting		✓	✓
	LAN Ethernet connection	AC	AC	 Animal weighing		✓	✓
	WiFi connection	AC	AC	 Dynamic difference weighing		✓	✓
	Bluetooth connection	AC	AC	 Free conversion Area Conversion		✓	✓
	External I/O (tare, print ..., relay)	OD	OD	 Statistics		✓	✓
	Integrated alibi memory (with EU type approval)	✓	✓	 Recorder		✓	✓
	Support IQ/OQ/PQ	✓	✓	 Minimum sample weight		✓	AC
	Integrated audit trail	✓		 Automatic repeatability test		✓	✓
	Below balance weighing	✓	✓	 Maintenance Service Package		✓	
	Anti-theft protection (mechanical and code)	✓	✓	 Net total application: Add up summation (Dosing program)		✓	✓
	Easy alphanumeric input	✓		 Checkweighing		✓	✓
	USB host (barcode, flash drive, ...)	✓	✓	 Density determination		✓	✓
	RS232 interface (for PC/printer)	✓	✓	 Buoyancy error suppression technology		AC	
	Direct tare and reference weight entry	✓		 Peak Hold		✓	
	Human Interface Device (for PC)	✓	✓				
	USB device (for PC)	✓	✓				
	Clock (GLP/GMP printout, auto-cal, etc.)	✓	✓				
	External Calibration, with calibration warning	✓	✓				
	Self-calibration system, fully automated	AC	AC				
	Identifier and label handling	✓					
	EU/OIML/UKCA type approval	AC	AC				

DIMENSIONS

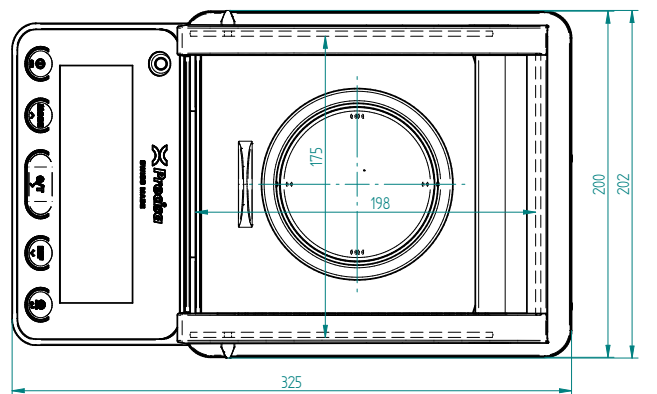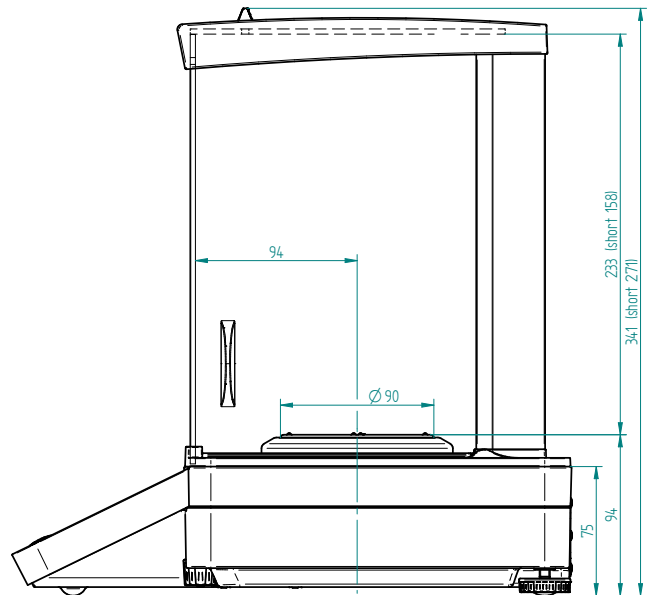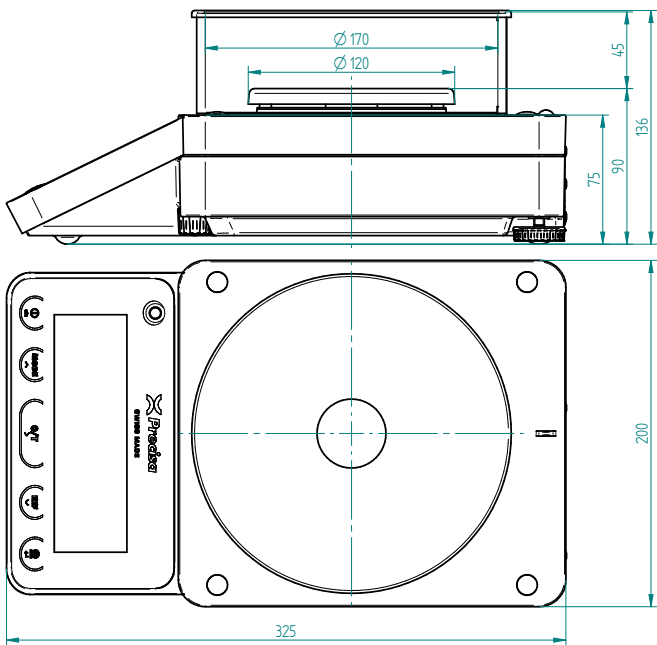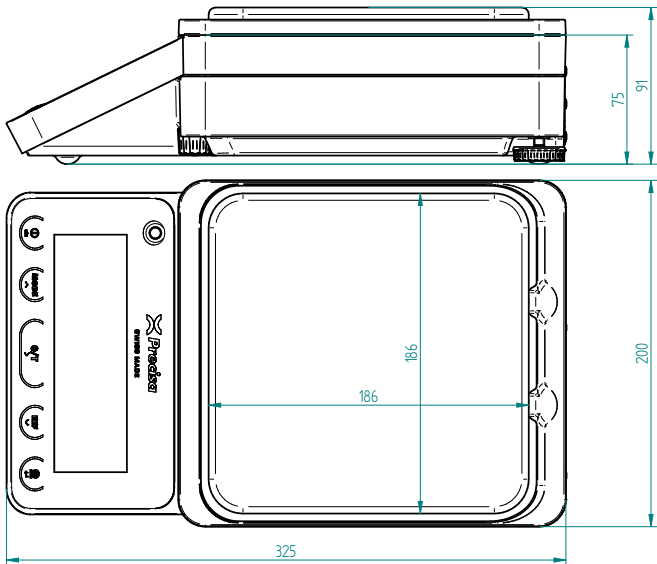

TECHCOMP GROUP

In addition to Precisa, Techcomp Europe comprises of the following companies:



Precisa Gravimetrics AG

Moosmattstrasse 32
8953 Dietikon, Switzerland
Tel. +41 44 744 28 28
Fax +41 44 744 28 38
Email: info@precisa.ch
www.precisa.com

ISO 9001
BUREAU VERITAS
Certification

N°CH4638370

