


# 350-8391-001

Bluetooth for Printer LP 4024,  
Power Supply via USB Host, 520

Bluetooth für Drucker LP 4024,  
Speisung über USB-Host, 520

## Scope of delivery – Lieferumfang

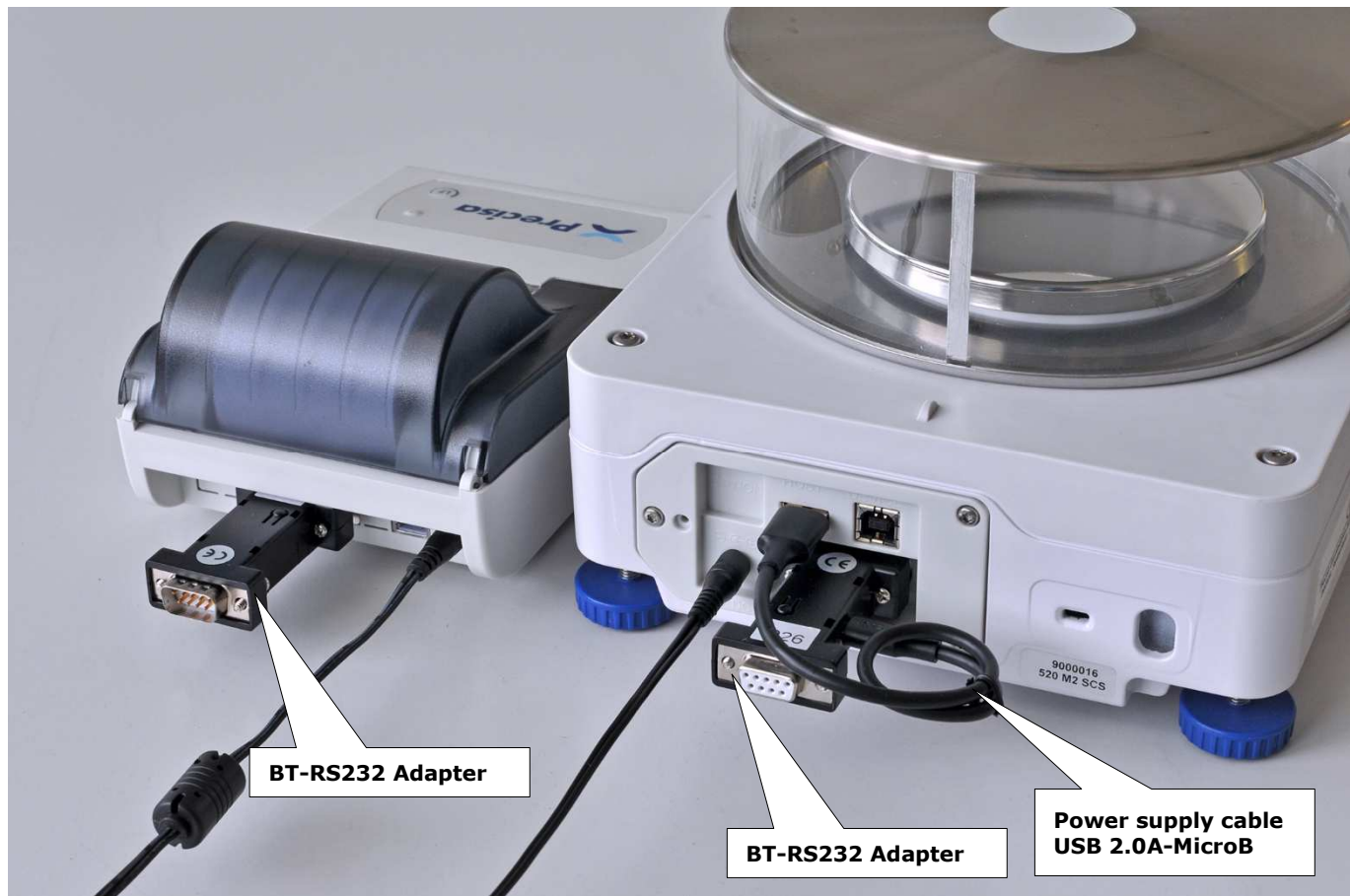

| Pos. | Pcs. | Item number  | Item                              | Artikel                        |
|------|------|--------------|-----------------------------------|--------------------------------|
| 1    | 1    | 350-8391-001 | BT-RS232 Adapter, pair            | BT-RS232-Adapter, Paar         |
|      | 1    |              | Pwr. supply cable USB 2.0A-MicroB | Speisungskabel USB 2.0A-MicroB |

## Starting up

The adapters are coupled ex works (9600Bd, 8-No-1; handshake is not operated).

1. Plug an adapter into the balance and power it from the balance's USB host using the included cable.
2. Plug the other adapter into the printer. This is powered directly from the printer.
3. Switch on balance and printer.
4. Ensure the serial interface of the balance is set to 9600Bd, 8-No-1.

## Inbetriebnahme

Die Adapter sind ab Werk gekoppelt (9600Bd, 8-No-1; Handshake wird nicht bedient).

1. Einen Adapter an die Waage stecken und mit beiliegendem Kabel über den USB-Host der Waage speisen.
2. Den anderen Adapter an den Drucker stecken. Dieser wird vom Drucker direkt gespeisen.
3. Waage und Drucker einschalten.
4. Sicherstellen, dass die serielle Schnittstelle der Waage auf 9600Bd und 8-No-1 eingestellt ist.