
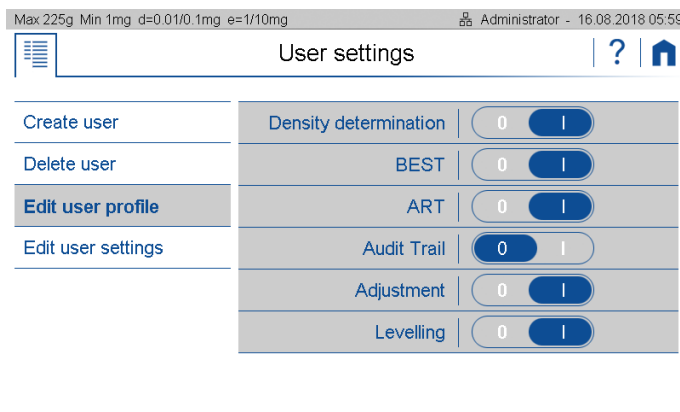
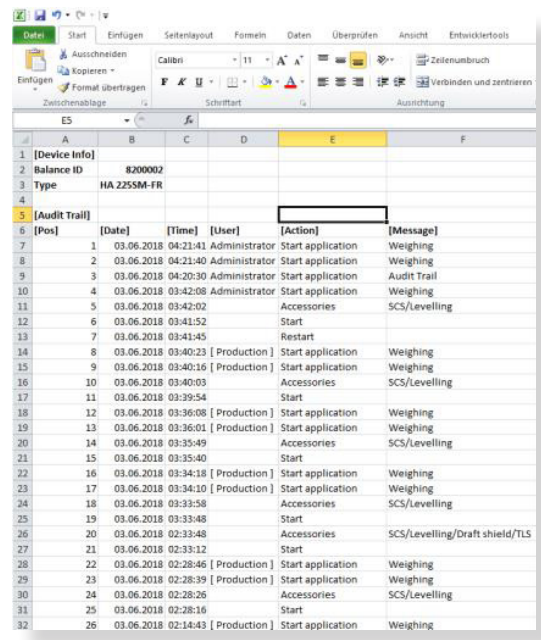


RELEASE NOTES SERIES 390

WHAT'S NEW IN VERSION 1.1.0


- Audit trail**

- Search function (by date)
- Manual entries possible
- Detailed display of the operations
- Export for detailed analyses on USB flash drive possible
- Access rights can be adapted individually via user management
- Free for all models

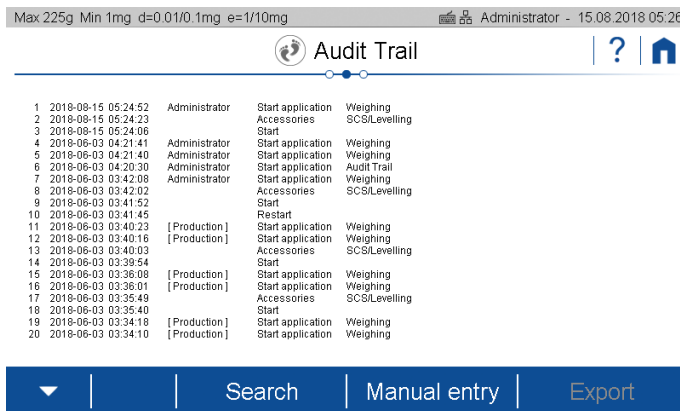


- **Density determination** extended

- Individually switchable air buoyancy and displacement correction
- Includes buoyancy due to density kit (Precisa or self-made) and air buoyancy

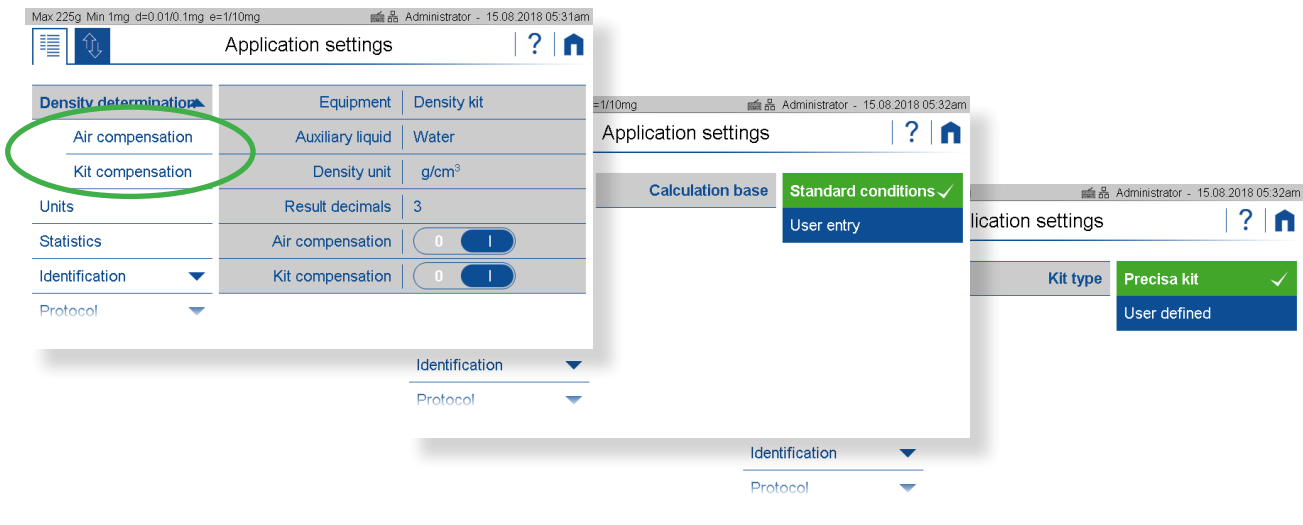


- In combination with the AOS, the error influence due to air buoyancy can be reduced to an absolute minimum (automatic measured value recording).
- Optionally, measured values from externally measuring instruments (pressure, temperature and relative humidity) can be entered manual.
- The workflow automatically requests the input of 3 parameters. The first time the workflow is used an input window is displayed. For the subsequent runs, the workflow points the possibility of manual input.

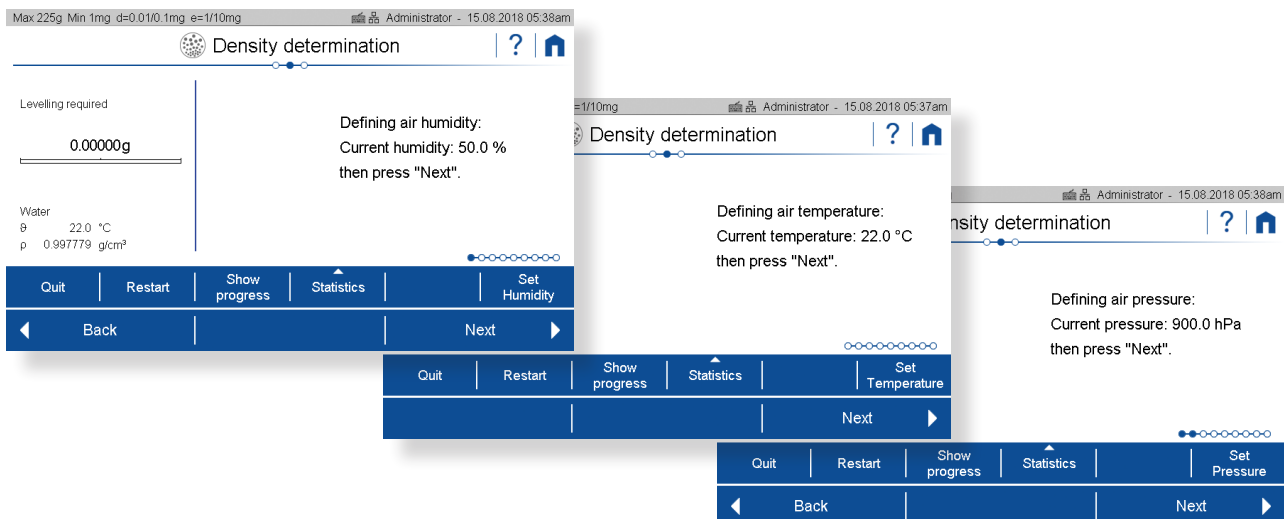


- **Statistic** extended

- Newly available in the density application
- In addition „ID’s“ can now be used

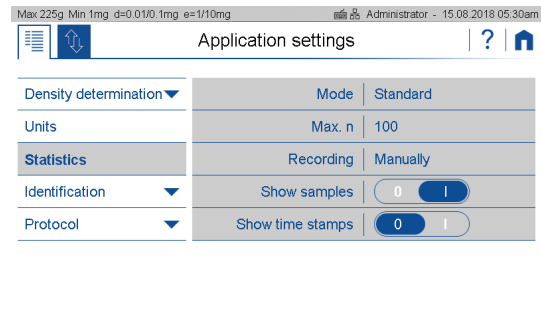


- **Auto-Standby**

- Screen saver with moving Precisa logo
- Including restart, abort selection

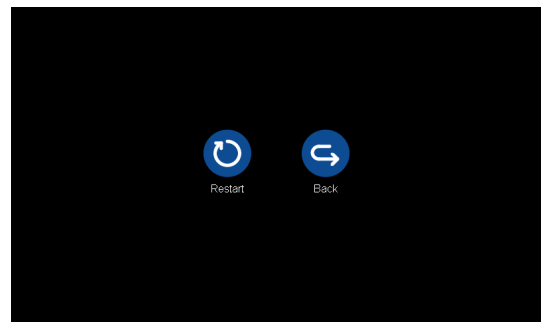


- **Statusbar**

- Icons for periphery added
- Icons are shown as soon as a connection is made (e.g. ethernet, mass storage, keyboard etc.)

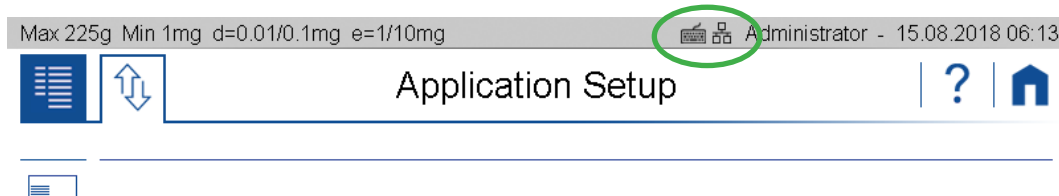


- **CAL button (Adjustment)**

- „CAL button“ selectable as soft button (customisation - via action bar in the active application)
- Adjustment can be started directly from the current application

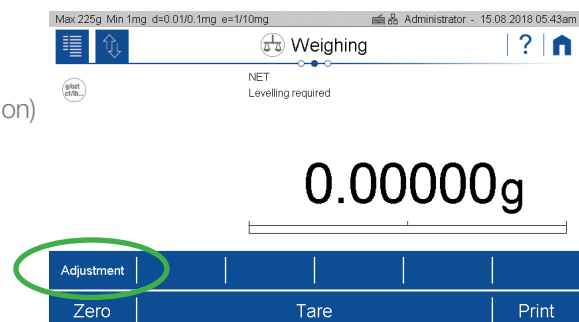


Folder soft buttons

- Actionbar can be used even better
- Easy to assemble via customizer
- Folder soft button can be named individually

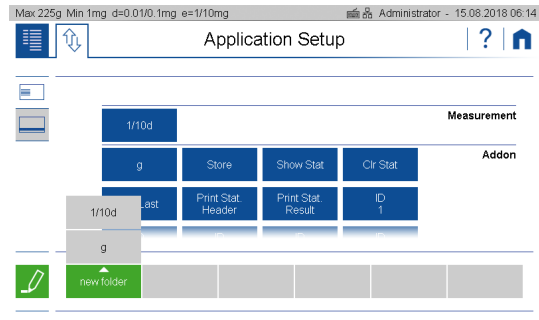

Date and time

- More efficient input mask to set the date and time


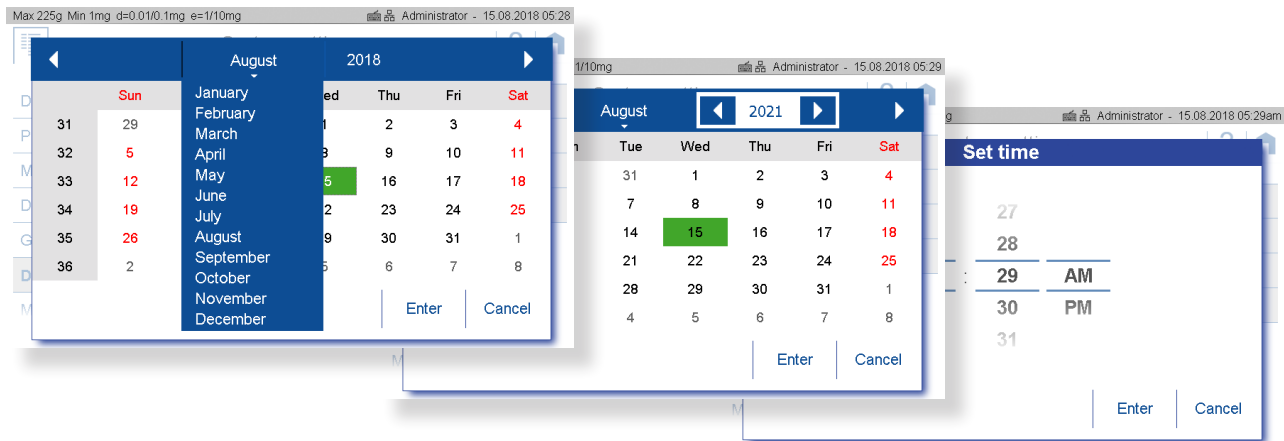


TLS

- Default settings for „Open draft shield“ new on „hold“


Units

- Fast unit switching with toggle button

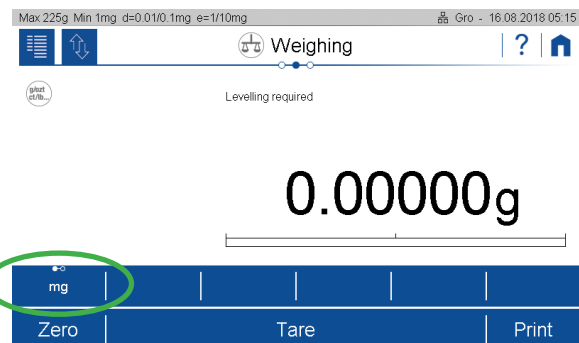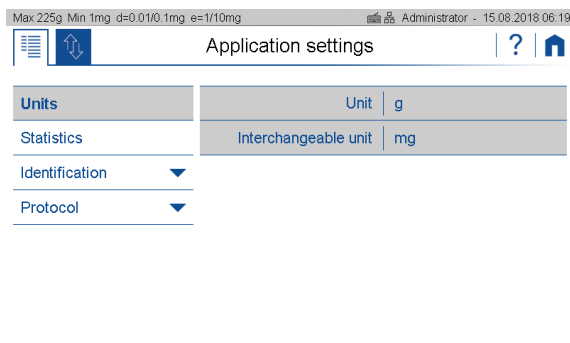

- **Logout**

- Fast logout implemented

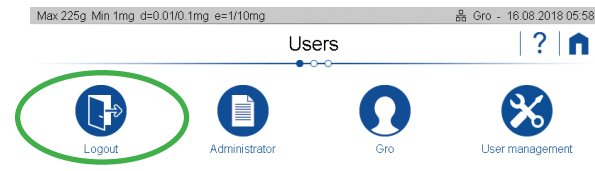

- **Software installation process** simplified

- Automated installation process for software update (silent installation)
- The installation process is initiated automatically by plugging in a USB flash drive

- **Minor bug fixes**